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**EN DIGITAL THERMOMETER model Thermo Base USER'S MANUAL**

- Measuring time: 60±10 seconds - oral, rectal; 100±20 seconds - axillary.
- Last measurement memory
- Fever alarm
- Auto shut-off
- Water proof


Thank you for choosing TM Gamma® Digital Thermometer Model Thermo Base. We are sure that having appraised worthy the high quality and reliability of this device you will become a regular user of the products of English Trademark Gamma®. Before starting to use the thermometer please study the instructions carefully.

**IMPORTANT SAFETY INSTRUCTIONS**

1. The thermometer is suitable exclusively for measurement of the body temperature!
2. The minimum measurement time until the signaling tone «beep» must be maintained without exception!
3. The thermometer contains small parts (battery, etc.) that can be swallowed by children. Therefore never leave the thermometer unattended to children.
4. Do not store the unit where it will be exposed to direct sunlight, dust or humidity. Avoid extreme temperatures, dropping or subjecting your thermometer to strong shocks.
5. Do not attempt to disassemble the unit, except to replace the battery.
6. Avoid ambient temperatures above 55 °C. NEVER boil the thermometer!
7. Store the unit in the protective case when not in use.

 **Batteries and electronic instruments must be disposed of in accordance with locally applicable regulations, nor with domestic waste.**

**FEATURES OF YOUR DIGITAL THERMOMETER**

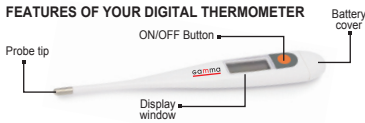**HOW TO TAKE YOUR TEMPERATURE**

1. Press ON/OFF button to activate. The unit would beep and display «188.8 °C». This LCD display test will go on for around 2 seconds.
2. Then, the temperature taken last time would be displayed with a «M» mark in the right corner.

3. Then, 37.0°C ("Lo" and "Co" flashing) will appear.
4. Put the thermometer to places you want to take temperature, into mouth, under armpit or inside the anus.
5. The beep will sound when temperature measuring is complete (when taking an oral and rectal temperature). If the signalling tone sounds and «C» is no longer flashing, this means that the measured increase in temperature is now less than 0.1 C in 16 sec. When taking an axillary temperature, ignore the beep sound and wait 3-5 minutes else before taking the reading.
6. Press ON/OFF button to turn off the unit. (The unit will automatically power off after approximately nine minutes after use.)

**MEASURING METHODS**

Measuring method	Position	Measuring time
<b>Taking an oral temperature</b> (in the mouth). Position the thermometer in one of the two pockets under the tongue, to the left or right of the root of the tongue. The measuring sensor must be in good contact with the tissue. Close your mouth and breathe evenly through the nose to prevent the measurement from being influenced by inhaled/exhaled air. Do not eat or drink cold or hot products before measurement.		60 ±10 seconds
<b>Taking a rectal temperature</b> (in the anus). This is the most reliable measuring method, and is especially appropriate for infants and small children. Apply a water-soluble lubricant to the probe cover. Gently insert the probe (MAXIMUM 1/2") into the rectum.		60 ±10 seconds
<b>Taking an axillary temperature</b> (under the arm). From a medical viewpoint, this method will always provide inaccurate readings, and should not be used if precise measurements are required.		100 ±20 seconds

**Note:** Rectal temperatures are generally 0.5°C /1.0°F higher than oral readings. Axillary temperatures are generally 0.5°C /1.0°F lower than oral readings.

**CLEANING INSTRUCTIONS**

Before and after each use, clean the thermometer with a soft cloth and ethyl alcohol diluted with water, or wash with COLD soapy water. Do not soak for an extended period of time or sterilize by boiling method, gas or steam autoclaves.

**REPLACING THE BATTERY**

The battery in your digital thermometer needs to be replaced when «□» appears in the lower right corner of the display. Replace the battery with a 1.55V, LR41 or SR41. To replace the battery:

1. Pull the battery cover away from the unit.
2. Use a non-metal pointed tool to remove the battery and insert the new battery with the positive(+)-side facing up.
3. Carefully slide the battery cover back into place.

**Note:** Please properly replace the battery away from small children and heat.

**SPECIFICATIONS**

Measurement range	32.0°C ~42.9°C
Low temperature display	Temperature <32.0°C: display: Lo°C
High temperature display	Temperature ≥42.9°C: display: Hi°C
Display resolution	0.1°C
Accuracy	34.0°C ~42.9°C ±0.1°C other range: ±0.2°C
Power Consumption	0.15 milliwatts in measurement mode
Normal temperature alarm (≤37.8°C)	Ten quick beeps for approximately 10 seconds
Fever alarm (>37.8°C)	Three quick beeps followed by a pause for approximately 10 seconds
Auto power off	Approx. 9 minutes
Battery	LR41 or SR41
Life of battery	More than 200 hours of continuous operation
Operating condition	Temperature: 10°C~40°C Humidity: 15%~95%RH
Storage condition	Temperature: -25°C~55°C Humidity: 15%~90%RH

\* The specifications and the outer design are subject to change for improvement reasons without previous notice.

**WARRANTY**

The digital thermometer is guaranteed for 3 years from date of purchase. The guarantee does not apply to damage caused by improper handling, accidents, not following the operating instructions or alterations made to the instrument by third parties. The guarantee is only valid upon presentation of the guarantee card filled out by the representative of the sales organization and its stamp and safety of the packaging.